


NOW ENHANCED!

Evoflex[®] Hip Belt


The most innovative hip belt on the market just got better! Our updated fastener kit gives you more control, and the new Evoflex Belt Extender improves versatility. What hasn't changed is how the Evoflex makes transfers easy with stiffened end straps that conveniently pivot out of the way when not in use. This versatile belt can be worn on the chest or legs as well as the hips. Users also love the secure feeling they get from the close-fitting, tapered straps.

Several buckle options to meet varied user needs, including the Rehab Latch[™] for those with limited hand function, the Small Push-Button for pediatrics, and the high-strength Swivel Buckle, preferred by active users.

Stiffened end straps rotate either direction and stay where you put them. Straps won't twist or fall into the wheels, avoiding dirt and damage. Convenient for sports or other intermittent uses as well as everyday support.

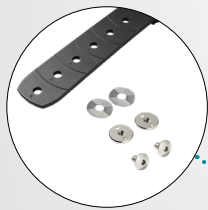

Low-profile fasteners fit into tight places. **NEW** wave spring washer allows variable tension to prevent unwanted rotation and help keep belt securely attached over time.

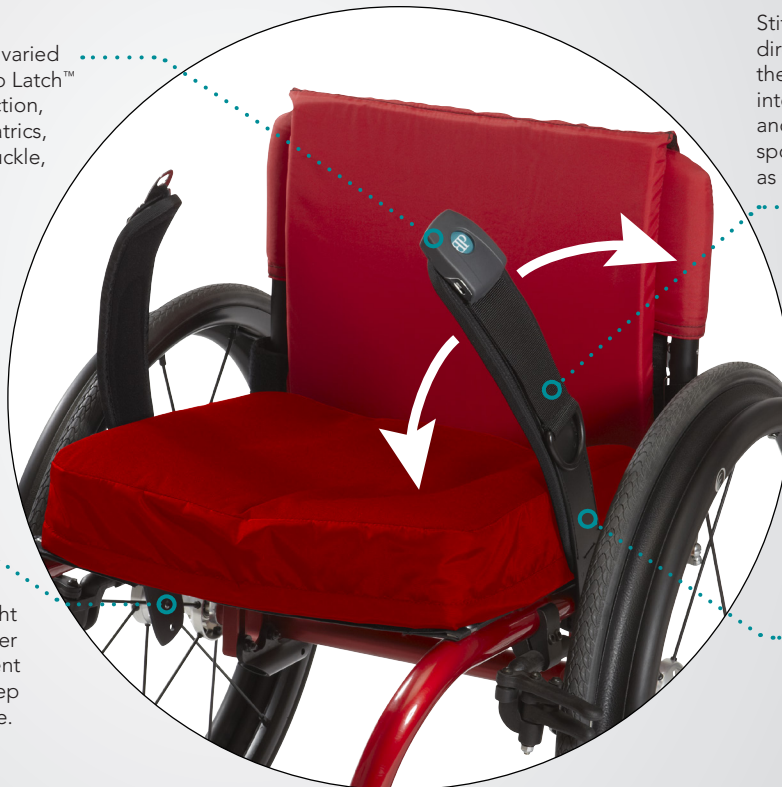

Patented, tapered strap design improves stability by providing a closer fit.

NEW! Evoflex[®] Belt Extender

Increase your Evoflex mounting range without compromising fit or performance. Designed specifically for the Evoflex Belt, our stiff, durable Evoflex Belt Extender lengthens the straps by up to 12.5 cm (5") on each side. Great for chairs with contoured seating, thicker seat cushions or distant mounting points. Buy separately or in a cost-saving bundle with the Evoflex Belt.



Made of high-impact plastic that maintains the stiffness of the Evoflex straps.

Smart locking pin adds a second connection point with the belt, preventing unwanted pivot. Alternatively, use the reverse side of the Belt Extender, without the locking pin, to allow pivot around obstacles.

Can be cut to desired length with a saw or bent to shape with heat as needed.




Visit us at bodypoint.com
206.405.4555 or 800.547.5716 | sales@bodypoint.com | Fax: 206.405.4556

© 2020 Bodypoint. All Rights Reserved. BMM328 2020.4



Measuring for the Evoflex®

To ensure proper fit and best performance, always select the Evoflex belt size that best fits the user. Then consider adding Evoflex Belt Extenders if longer straps are needed to reach the desired mounting point.

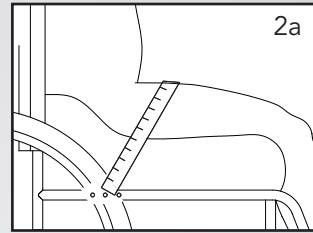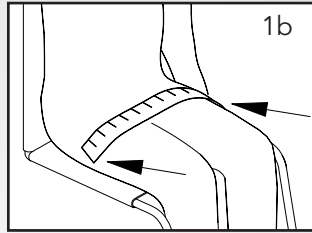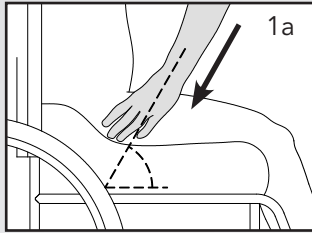


1. Determine the correct belt size for the user:


- Based on clinical goals, use arms to place user in desired position and determine the most effective position and angle for the belt. (Fig. 1a)
- At the selected belt angle, measure the distance over the body with a flexible tape to find the optimal length for the belt's padded section. (Fig. 1b)
- In the sizing table, refer to Column B to find the appropriate belt size based on padded length. If user is between sizes, consider factors such as growth, weight variations and clothing changes. However, for proper fit and function, pads should not extend past the user's body.

2. Determine total length needed for mounting:

- At the selected belt angle, measure from the mounting point on one side across the user's body to the mounting point on the other side. (Fig. 2a)
- In the sizing table, refer to Column D to find the maximum mounting length for the selected size belt.
- If belt will not be long enough to reach the mounting points, Evoflex Belt Extenders can be used to add up to 12.5 cm (5") length on each side. Both the Evoflex Belt and Extender can be cut to length if needed. Belt Extender can also be heat bent.

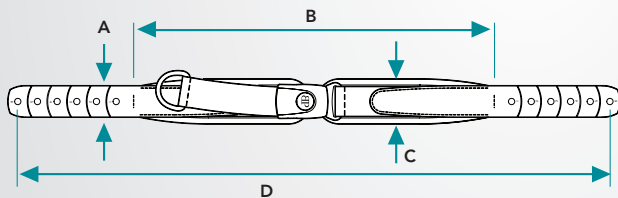

Size	A	B	C	D
S38	2.5cm (1")	38cm (15")	5cm (2")	63cm (24-3/4")
M46	4cm (1-1/2")	46cm (18")	5.5cm (2-1/4")	76cm (29-3/4")
L62	5cm (2")	62cm (24-1/2")	7cm (2-3/4")	96cm (37-3/4")

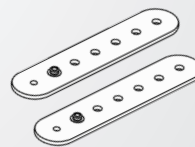

Belt Size	Max mounting length with Belt Extender
S38	88cm (34-3/4")
M46	101cm (39-3/4")
L62	121 cm (47-3/4")

Part Numbers

Evoflex®

Size	Push-Button	Rehab Latch™	Swivel Buckle	Small Push-Button
S38	-	-	EB275-S38	EB235-S38
M46	EB205-M46	EB225-M46	EB275-M46	-
L62	EB205-L62	EB225-L62	-	-

NEW! Evoflex® Belt Extender

One size is compatible with all Evoflex belts. Adds up to 12.5 cm (5") length per side, with holes in 25 mm (1") increments. Includes fasteners. Sold in pairs.



Buy together and save!

Evoflex® + NEW! Belt Extender

Size	Push-Button	Rehab Latch™	Swivel Buckle	Small Push-Button
S38	-	-	EB275-S38-EBE	EB235-S38-EBE
M46	EB205-M46-EBE	EB225-M46-EBE	EB275-M46-EBE	-
L62	EB205-L62-EBE	EB225-L62-EBE	-	-

Evoflex Hip Belts and Belt Extenders include low-profile fasteners for mounting; no additional end-fittings are required. Combine with Bodypoint's other modular belt mounting hardware as needed to meet installation needs. Not compatible with the FrameSaver clamp.

Visit us at bodypoint.com
206.405.4555 or 800.547.5716 | sales@bodypoint.com | Fax: 206.405.4556

© 2020 Bodypoint. All Rights Reserved. BMM328 2020.4

