

Healthcare
furniture,
enhancing
lives.

Caremed **Alrick.**

CSS Concave Mattress

Helps Prevent Falls And
Provides Safety



This is a deluxe grade mattress that features concave sides, so the patient is encouraged to lie in the middle of it, thereby minimising the possibility of falling. It facilitates pressure distribution, and is vapour permeable. The heavy-duty base reduces the chance of wear due to sliding. It is available in single and king single widths.

Specifications

Overall dimensions	Single	1960 x 860 x 150mm
	King Single	1960 x 1050 x 150mm
Product codes	Single	CSS
	King Single	CSSKS
SWL	200kg	
Warranty	2 yr	


Features

CORE

- ☒ Deluxe Grade Core Base with Deluxe Grade Inbuilt Convolute Overlay
- ☒ Deluxe Grade Concave Sides encourage patient to lie in centre of mattress and help prevent patient falling off mattress

COVER

- ☒ Combination Cover
- ☒ Upper 2 Way Stretch Assists with Pressure Care and Pressure Distribution
- ☒ Vapour Permeable Upper to draw moisture away from skin
- ☒ Heavy Duty PVC base provides durable base to promote cover longevity and minimize wear caused by mattress sliding on bed base
- ☒ Flap over zip helps prevent fluid penetration

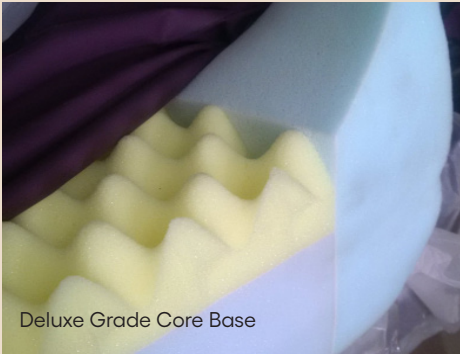

Deluxe Grade Core Base



Vapour Permeable Upper