


Healthcare
furniture,
enhancing
lives.

Caremed Alrick.

M9500 Bariatric 500kg
Mattress


Vapour permeable upper to
draw moisture away from skin
& impressive SWL of 500kg

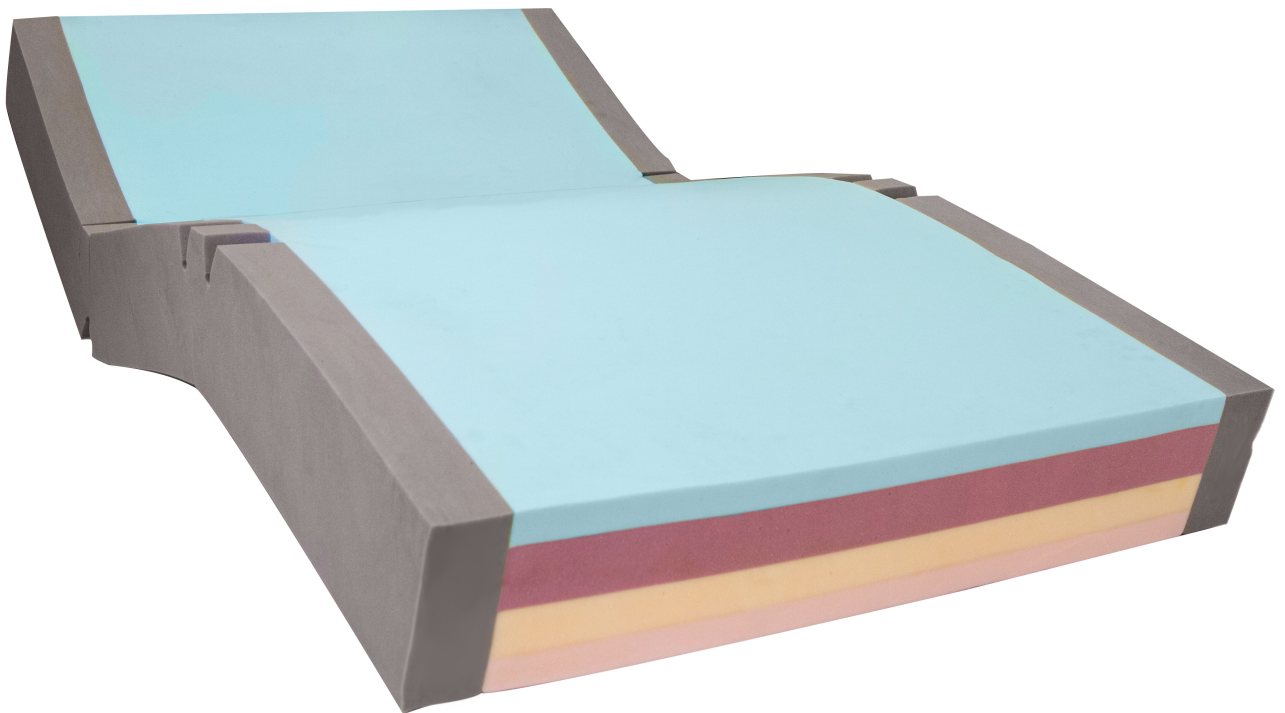

The M9500 Bariatric Mattress provides care, pressure relief and comfort for bariatric patients weighing up to 500kg. Through the use of proven foam combinations, unique construction methods and high performance cover materials, this mattress is the fitting choice for clinicians and facilities looking to achieve comfort, longevity and effective patient care.

Specifications

Overall Dimensions	King Single	1960 x 1050 x 175mm
	Extra Wide	1960 x 1150 x 175mm
	Double	1960 x 1350 x 175mm
	Queen	1960 x 1480 x 175mm
Product codes	M9500KS	King Single
	M9500EW	Extra Wide
	M9500D	Double
	M9500Q	Queen
SWL	500kg	
Warranty	2 yr	

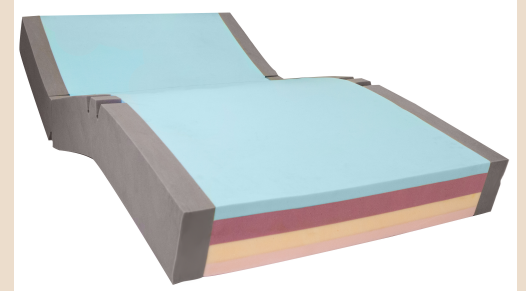


✓ 4-LAYERED GRADED PRESSURE CARE CORE

- High specification foam for longevity and performance
- Provides graded postural support and patient comfort
- Foam conforms to AS/NZS 4088.1 (Int):1993 and FMVSS302 & SAE T2002
- Anti-microbial treatment applied to foam maximises longevity and core micro-climate



✓ PREMIUM VISCO-GRADE, UPPER IMMERSION LAYER

- Visco-elastic cell structure provides low-shear surface, effective immersion and interface pressure relief
- Provides exceptional comfort
- Delivers postural support for patient repositioning and transfer

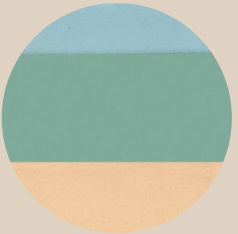

✓ PREMIUM POSTURAL SUPPORT SECONDARY LAYER

- Graded support and patient comfort in supine, prone, lateral and fowler positions
- Excellent core stability to encourage centralised patient positioning to minimise falls risk, entrapment and enhance pressure relief



✓ LAYERED FIRM BASE

- Provides protection against 'bottoming-out' during transfer, repositioning and whilst in the inclined/fowler position
- Graded for ultimate patient support



✓ TRIPLE SECTION, SPLIT BASE

- Enhances profiling with adjustable bed surface positions minimising excess strain on bed mechanisms

