Healthcare furniture, enhancing lives.

Caremed Alrick.

Total FlowPlus





CM-M-OLB Total Flow Plus – a replacement alternating pressure relieving mattress system with 5" air cells and a self-inflating 2" sealed foam base forming an in-use height of 7" for effective prevention and treatment of patients at high risk of developing a pressure ulcer in hospital, nursing and care home environments.

With simple to use settings and functionality, the Total Flow Plus enables the care provider to quickly set up the mattress system and have it operating at optimal pressure for a specific user with ease.

The wide range of features make the Total Flow Plus a highly versatile and cost effective solution which has been manufactured to comply with the most stringent quality and in-use guidelines.

Total Flow Plus incorporates cutting-edge Total Ait Technology - constant auto pressure sensing through an air sensor and microcomputer regulator to provide ultra-low interface pressure and promote an even ambience for natural sleep patterns and stable blood circulation with silent pump operation.

TECHNICAL SPECIFICATION

Mattress Specification

- Mattress Dimensions: Single 200 x 90 x 12.7 cm + 5cm foam base King Single 200 x 107 x 12.7 cm + 5cm foam base
- > Mattress Weight: 5.5kg
- > Cell Structure: Cell-on-cell
- > CPR: Quick release toggle
- > Max Load: 200kg
- > Heel function (3 x heel cells with one way valves)

Pump Specification

- > Operating System: Total Ait Technology
- > Pump Dimensions: 34.5 x 24 x 12.5 cm
- > Pump Weight: 2.8kg
- > Cycle Time: 10, 15, 20, 25 minutes
- > Memory function
- > Air Output: 8 litres/minute
- > Power Input: AC220-240V, 50-60Hz
- > Power Consumption: 7W
- > Fitted with hanging brackets

Cover Specification

- > Multi-stretch PU
- > Vapour permeable
- > Water resistant
- > Machine washable
- > Antimicrobial
- > 360 degree zipper
- > White underside of cover



Nylon TPU coated cells



Independently removable cells



Covert sealed air cells



Removable heel cells



Pillow function




CPR release



Anchoring straps



Zipped storage bag

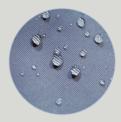

Multi-stretch, waterproof & vapour permeable PU cover



Micro-perforated cells



Antimicrobial. machine washable



360 Degree zipper



Robust pump connection

Silent pump operation

to avoid disturbance



Transport mode

Multiple cycle

time control


modes

CONTROL PANEL

Auto pressure sensing for consistent distribution of air

Visual and audible warning of air loss. power failure and alternate failure

Control panel lock to prevent unauthorised Or accidental tampering (hold for 2 sec)

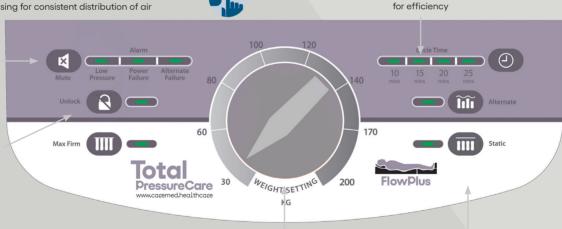

Weight setting Static/alternating control knob

Clinical pressure mapping evaluation

Australia

VIC 2/77 Ricky Way Epping VIC 3076

NSW 104 Loughnan Street Lake Cargelligo NSW 2672

T 1800 282 001 sales@caremed.healthcare www.caremed.helthcare

United Kingdom

Port Stanley House Falkland Close Coventry, CV4 8AU

E & O.E.

T 01926 889677 updates@caremed.healthcare www.caremed.healthcare