

# Unique solution for safe patient showering, bathing and transfer

The versatile Pacific Shower-Bathing Trolley provides a safe environment for showering, bathing and transfer. The soft, padded top is ideal for patients with delicate skin or disabilities.


The electric high-low and unique electric tubsides actuation provides an ergonomic working height for caregivers and ensures safer patient transfers.

The Pacific Shower-Bathing Trolley has been carefully designed to be intuitive and easy to use meaning minimal user training is required. The one piece seamless design also makes it easy to clean.

# **FEATURES**

- Fully electric controls with battery backup
- Tub sides lower to become flush with patient platform
- High performance LINAK and TENTE components
- 200kg Safe Working Load
- 125mm corrosion resistant castors
- Safe 24V DC system with remote battery charger

- Power cords available to suit different countries
- Soft, padded polyurethane tub
- Low battery warning light on handset
- Waterproof pillow
- Push handles at both ends
- Siderail (option)

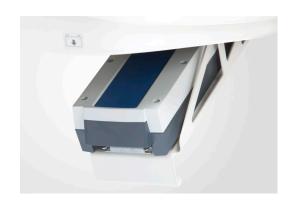









## Doorway width (mm)

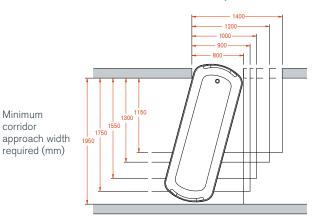

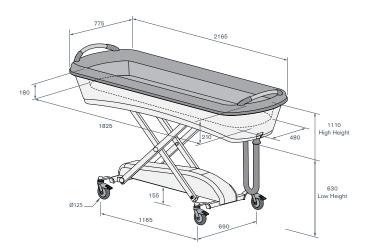

## **SPECIFICATIONS**

Minimum corridor

required (mm)

MODEL:	ST 2.1 and ST 2.2 (Base raised 30mm)	
IP RATING:	IPX5	
EEC CLASSIFICATION:	Class 1, Non-invasive	
MAXIMUM SWL:	200kg (incl. patient, accessories & water)	
TROLLEY MASS:	100kg	

- · Unique, seamless deep tub design enables immersion bathing
- Fully electric controls enable a comfortable working height adjustment
- · Easy cleaning for improved infection control
- · Safe 24V DC system with remote wall mounted battery charger
- Ergonomic push handles at both ends
- Castors: 125mm (synthetic with sealed wheel bearings), 3 individual locking, 1 directional locking
- 360° bumper
- · LINAK controls with 2 actuators, control box, handset & battery

# **ELECTRICAL SPECIFICATIONS:**

- Power source: Rechargeable battery, 24V DC, 2.9Ah
- Battery charger input: 110 or 240V AC, 50-60 Hz, 650mA (max)
- · Battery charge time: 6 hours
- Number of full load cycles: approx. 30
- Time to raise or lower patient platform: 30 sec
- Time to raise or lower tub sides: 22 sec

# SAFETY COMPLIANCE

• Developed in accordance with the requirements of AS/NZS 3200.2.38:2007 & IEC60601-2-38:1996/10

## **DIMENSIONS**

- Trolley overall: L: 2165mm, W: 775mm • Tub platform L: 1825mm, W: 480mm • Wheel base: 1165mm, wheel track: 690mm Castor diameter: 125mm single wheel
- Maximum tub depth (drain end): 210mm
- Maximum tub depth (head end): 180mm

ADJUSTMENT RANGE:	ST2.1	ST2.2
Maximum overall height (sides up, incl handles):	1155mm	1185mm
Maximum trolley height (sides up, excl handles):	1080mm	1110mm
• Minimum overall height (sides down, incl handles):	600mm	630mm
• Minimum trolley height (sides down, excl handles):	525mm	555mm
Maximum patient platform height (drain end):	870mm	900mm
Maximum patient platform height (head end):	900mm	930mm

## UNDERSIDE CLEARANCES (FOR PATIENT LIFTERS)

• Maximum floor clearance: 125n	nm 155mm
---------------------------------	----------

# www.howardwrightcares.com

Paraite Road, Bell Block PO Box 2786 PO Box 3003, Fitzroy Taren Point New Plymouth 4341 New South Wales 2229 New Zealand Australia T. +64 (6) 755 0976 T. 1800 120 727 F. +64 (6) 755 0908 F. 1800 120 717



Design and Development, Production, Sales, Distribution and Servicing of Medical Beds, Stretchers and Accessories

Howard Wright Limited has a policy of continuous product improvement and reserves the right to change product specifications without notice. Howard Wright Limited is certified to ISO13485:2012.





