



Little Wave "Flip"

US Pat# 9,801,766

FLIP THE SCRIPT ON ORDINARY.

Designed to be as unique as your child.

Flip is designed for the smallest of riders with some of the biggest advancements in the industry! From its flat folding design to its use of TAPER LOK Technology the Flip is simply the most advanced pediatric tilt-in-space available today.


Looks like a stroller.
Works like a wheelchair.

The Little Wave Flip integrates design elements from strollers into the durability of a wheelchair.



Express your personality!

With over 500 different color combinations, colored mag wheels, and designer spoke guards; it's easy to personalize your Flip.

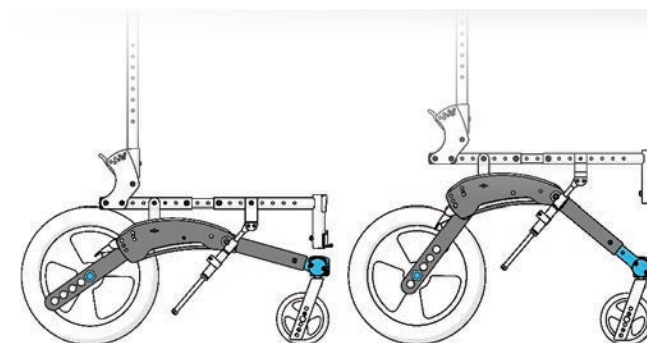

Don't struggle with growth.
Embrace it.

The Flip offers a new paradigm in adjustability. With fewer parts than any other chair, seat height, depth and width have never been so easy to change.



Attendant Wheel Lock

A hands-free, easy to use stroller like wheel lock option that makes a parent's life just a little simpler.



Flip Up Anti-Tips

Tired of struggling with anti-tips? Our design allows the anti-tips to stay connected while taking on every day obstacles, such as curb climbing.




Drum Brake

We are not eliminating the wheel lock, instead we are adding the unique option of a true brake.



Tough, Strong & Easy to Fold

With the Flip of a lever you are able to activate the **Ergo Folding Back** for a controlled, effortless fold.

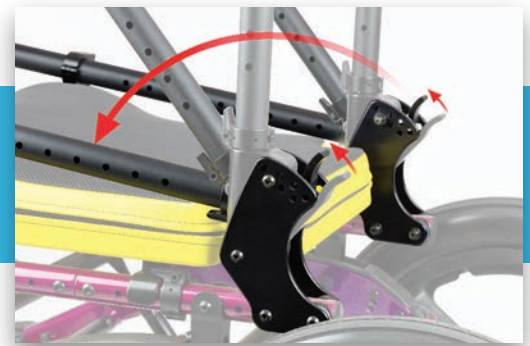

SPECIFICATIONS

Transport Weight:	21 lbs.*	9.53 kg*
Seat Width:	10-18 in.	25-46 cm
Seat Depth:	12-20 in.	30-51 cm
Seat Height:	13.5-20 in.	34-51 cm
Tubing Size:	1 in.	2.54 cm
Weight Capacity:	165 lbs.	75 kg
WC-19 Transit:	Yes	Yes

** 14"x14" Frame with Attendant Wheel Locks.
Armrests, Rear Wheels and Footrests Removed.*