



ZOOM!
around your world

The Zippie ZM-310™ combines our smallest, most reliable power wheelchair base with clinically proven pediatric seating and positioning options to provide children maximum accessibility and independence.



ZM-310



COMPACT



TILT & LIFT



POWER



ZM-310

OUR SMALLEST

Mid Wheel Drive Base

The Zippie ZM-310's compact base is only 24" wide and 35" long for accessing tight indoor spaces. Its mid-wheel-drive design is intuitive to drive and turns on a dime with a turning radius of only 19". The combination allows your child to drive effortlessly through narrow hallways, between desks, and beside friends and family.

13" DRIVE WHEELS

The ZM-310 features 13" drive wheels with solid tires for climbing obstacles and gliding over uneven terrain while remaining secure, comfortable, and flat-free.

POWERFUL

Motors & Batteries

Our reliable and high performance motor package features 4-pole, 10kph motors that keep up with your child's most playful and energetic days. The ZM-310's long-lasting batteries provide the range needed to travel longer distances on a single charge.



ELECTRONICS

Specialty Input Devices

When a standard joystick is not the best option for your child, Specialty Input Devices help make independent mobility possible. Quickie and specialty input device leaders Switch-It® and ASL® offer devices such as switch controls, head array, chin controls, and sip and puff that focus on the strengths and capabilities of each child to create a custom driving platform they can operate with ease.

PG Drives VR2 Electronics

PG Drives standard, non-expandable VR2 electronics are easy to use, offer simple control of common functions, and use a toggle switch to control power tilt.

PG Drives RNET Electronics

PG Drives expandable R-net electronics offer more advanced programming including multiple drive profiles, which allow speed, acceleration settings to be preset for different environments such as school, home, and R-net is compatible with multiple joysticks, displays, and specialty input devices.





POWER TILT

50° of power tilt allows you or your child to change their position and shift weight from the pelvis to the back to reduce the risk of skin breakdown. For children with limited pelvic, thoracic, or head control, this change in position helps decrease the potential of developing postural deviations down the road.

The ZM-310 also offers **50° of manual tilt** including a folding handle for easy operation.



POWER SEAT ELEVATE

9" or 12" of power seat elevate helps your child independently access shelves, cabinets, and light switches and provides the opportunity to interact with their peers at eye level. Independence at last!



MANUAL RECLINE

40° of manual recline also allows you to change your child's position or posture to assist with feeding, digestion, respiratory function, and visual orientation.

PEDIATRIC SEAT SIZES

The ZM-310's ASAP II (Adjustable Seating & Positioning) seat frame has been designed specifically for pediatrics and is available in sizes ranging from 10" x 10" to 16" x 19", of course, it adjusts to accommodate growth.



ANGLE ADJUSTABLE CENTREMOUNT FOOTREST

Angle Adjustable Centremount Footrest with Pediatric Footplate offers 3" to 13" lower leg lengths for the best fit and close access.



ZM-310 ALLOWS
YOUR CHILD
TO NAVIGATE
TIGHT SPACES
AND IS CAPABLE
OF TACKLING
TOUGH TERRAIN

JAY®



JAY® SEATING OPTIONS

The JAY Zip™ Back, JAY J3™ Back, any JAY cushion, and various JAY Positioning Supports may be ordered on-chair with the ZM-310 and JAY custom seating systems such as SureFit™ are available through Compatibility Assurance.

This guarantees compatibility and prevents delays caused by incompatible components.

ACCESSORIES

Complement your Zippie ZM-310



Fixed Vent Tray



Active Controls
Power Buddy



RAM® Phone
Mount



RAM® Tablet
Mount



RAM® Cup Holder



Active Controls
ReCharge™



Switch-It® Remote
Stop Switch



O2 Holder



Cane/Crutch Holder




Walker Holder

The ZM-310 is available with 14 shroud color options!


IMPORTANT NOTE: Print processes cannot accurately represent paint finishes. Actual colors may vary.



Apple Green

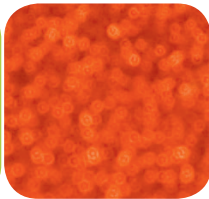


Copper

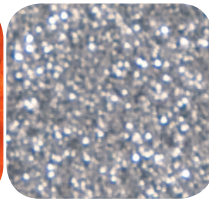

Silver




Matte Black

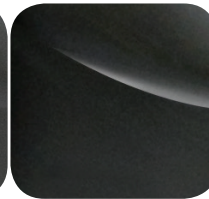


High Gloss Black



Yellow

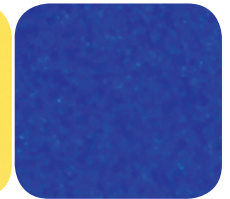

Candy Blue

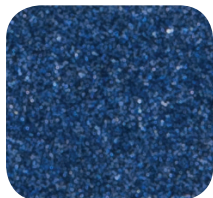

Midnight Blue



Carbon Fiber



Pink Camo



Pearl Pink



Candy Apple Red



Pearl White



Grunge Skulls

ZIPPIE ZM-310 SPECIFICATIONS

Drive Wheel	Mid Wheel	Drive Wheels	13" with flat-free insert
Max Speed	10 kph	Batteries	50 AmpH
User Weight Capacity	75 kgs	Tilt	50° Optional Power or Manual Tilt
Electronics	PG Drives VR2 or R-Net	Recline	40° Optional Manual Recline
Width	10" to 16"	Seat Elevate (Lift)	Optional: 9" to 12"
Depth	10" to 19"	Vent Tray	Fixed or Lift-off
STFH	15.5" to 21.8"	Transit Tested	WC - 19
Overall Base Width	24"	Warranty	Frame - 5 Years / Parts - 1 Year
Overall Length	35"	TGA	100375
Turning Radius	19"		