



CAT I™
standing frame



Standing frame **CAT I™** was created with a view to help children in adoption of the upright position. The upright standing contributes to the acquisition of a new experience by children from the nearest environment and stimulates their psycho-physical development. In addition, it regulates the physiological processes and oxygenates the body of the child. There are many advantages during using of the **CAT I™**

Adjustable tilt angle (from **25°** to **85°**) even with a child in the device. Extensive system of belts and supports makes the child feel safe and secure.

As far as this product is concerned, it is made of ecologic and environmentally friendly materials, and is also hand-painted. **CAT I™** is available in 1 size.

IT IS RECOMMENDED IN THE FOLLOWING CASES:

- Cerebral palsy – various types and forms.
- Brain and spinal cord injuries with lower limbs paralysis and weakness, balance and motor coordination disorders with different etiology.
- Demyelinating disease (e.g. multiple sclerosis – SM).
- Meningocele.
- Muscles diseases with weakness and paralysis (dystrophy, myasthenia gravis, myositis).
- Genetic disorders.
- Degeneration and neurological diseases.
- Other conditions connected with muscles weakness, paralysis and musculoskeletal disorders.
- After head and spinal cord injuries and spine trauma with paralysis and paresis.



STANDARD EQUIPMENT:

Frame with castors and brakes | Foot plate with foot support and belts | Gas spring | Tilt-in-space | Chest and pelvic pads with safety belts | Chest and pelvic supports | Tray with plate | Tray pad

Our passion – your independence!

Colors may differ from the presented images. The manufacturer reserves the right to make changes in products without updating the directory.
Current offer on our website.

Akces-MED sp. z o.o.
Jasionka 955B
36-002 JASIONKA

export@akces-med.com
+48 (17) 864 04 77
www.akces-med.com

Version: 14/05/2018



ADVANTAGES:

**Stable structure**

The product ensures adequate stability and safety of a Patient.

**Eco-friendly materials**

We use only human and environment friendly materials.

**Quiet operation**

The device is equipped with breathable, hypoallergenic 3D spacer fabric. Now you can be sure the skin stays dry for longer!

**Breathable upholstery**

The device is equipped with breathable, hypoallergenic 3D spacer fabric. Now you can be sure the skin stays dry for longer!

**Easy to use**

Intuitive way of adjustment and easy assembly.

**Pneumatic support**

Gas spring helps in proper tilt angle adjustment of the device.

**Footrest with 3-dimensional adjustment**

Footrests with 3D adjustment allows for smooth foot setting in all directions.

**24 months warranty**

We offer 24 months warranty on the product. We also provide after sales service.



Gas spring helps you to set the proper tilt angle of the device.



System of belts and supports makes your child feel safe and secure.



Tray with storage for toys.





BREATHABLE UPHOLSTERY ATMO™



BREATHABLE



EASY TO CLEAN



DURABLE



ECOLOGICAL



ANTIALLERGIC

ANTIMICROBIAL
& ANTIFUNGALGrass
StandardSea
Optional*Latte
Optional*Red
Optional*

WATERPROOF UPHOLSTERY HYDRO™



WATERPROOF



EASY TO CLEAN




DURABLE



ECOLOGICAL



FIRE RETARDANT

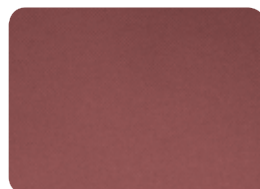Green Grass
Optional*Blue Sea
Optional*Brown Sugar
Optional*Hot Red
Optional*

MEDICAL UPHOLSTERY EXTREME-MED™

BPA FREE
PHTHALATE FREE300,000
MARTINDALE
RUBS
ABRASION

DURABLE




ECOLOGICAL



CRIB7

ANTIMICROBIAL
& ANTIFUNGAL

STAIN RESISTANT

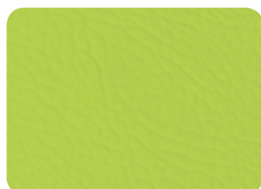Citrus Green
Optional*Candy
Optional*Pillarbox
Optional*

MEDICAL UPHOLSTERY MEDICUS™



WATERPROOF



EASY TO CLEAN



STEAM PERMEABLE



ECOLOGICAL





FIRE RETARDANT

MED Sea
Optional*MED Blue
Optional*MED Skin
Optional*

Optional* colors available for special request. Lead time ca. 14 days.

*Our passion - your independence!*Colors may differ from the presented images. The manufacturer reserves the right to make changes in products without updating the directory.
Current offer on our website.Akces-MED sp. z o.o.
Jasionka 955B
36-002 JASIONKAexport@akces-med.com
+48 (17) 864 04 77
www.akces-med.com

Version: 14/05/2018

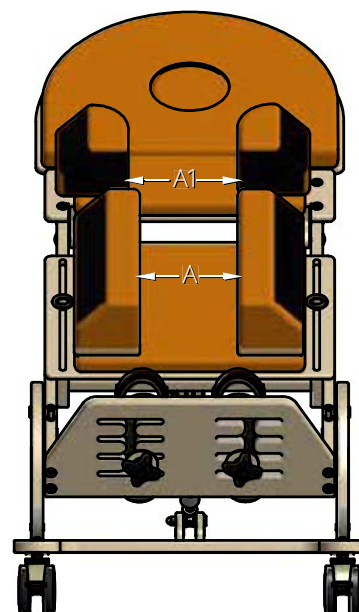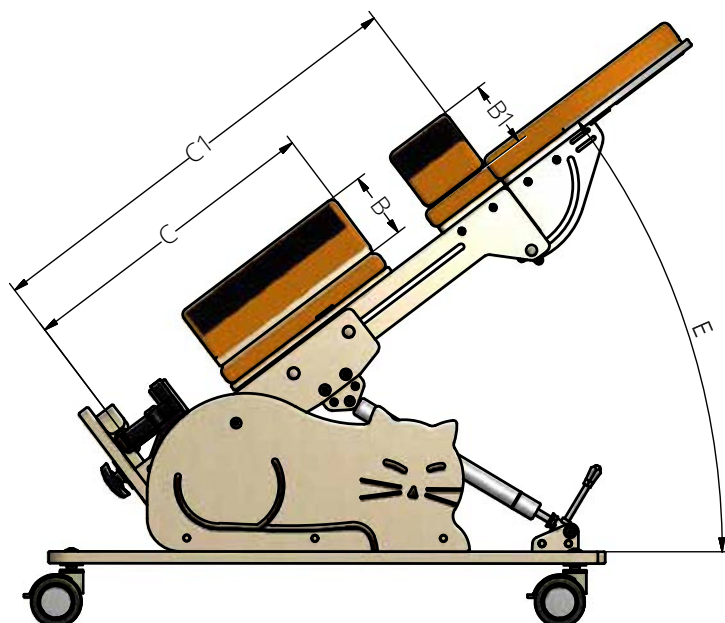

Symbol	Measurement	Unit	Size 1
A	Pelvic supports width	[cm]	12-23
B	Pelvic supports depth	[cm]	12
C	Pelvic supports height adjustment	[cm]	38-60
A1	Chest supports width	[cm]	12-21
B1	Chest supports depth	[cm]	12
C1	Chest supports height adjustment	[cm]	50-72
	Footer length	[cm]	17
E	Tilt angle	[cm]	25-85
	Max. user weight	[cm]	25
	Indicative user height*	[cm]	70-100

Product dimensions & weight		Unit	Size 1
Width		[cm]	54
Length		[cm]	100
Height		[cm]	73
Weight		[kg]	17

Note: The value cannot be main factor for selection of product size.
Note: Measurement tolerance range is +/- 1 cm / ½ in.

**Safety****TANGO Wheel (75/100/125mm)**

Noiseless, easy to maintain, does not leave marks, does not corrode.

**Stabilization****3D foot adjustment**

They allow for a smooth foot setting in all directions thanks to the use of the hinge ball. You can adjust your foot setting on the foot plate (front/back, left/right) and also set your foot in space (plantar/dorsal flexion, pronation/supination of the foot).

**Head supporting belt**

Attachable to the headrest with Velcro straps. Indicated for people with muscle weakness to prevent falling of the head.

**Comfort****Tray for manual therapy**

The tray allows active exercise assisting upper limbs with a particular focus on a palm function.

**Tray with bowl**

More toys or other therapeutic items can be put into this bowl now.