



**KIDOO™**  
positioning chair



**KIDOO™** is a therapeutic chair that is perfect for therapy and fun as well as other daily activities such as learning or having meals. The chair can be used at home, school, kindergarten or in a group of peers. **KIDOO™** turns out to be an indispensable tool in the rehabilitation process, both in terms of therapy and social integration.

This therapeutic chair provides the child with proper support and safety, moreover, it can be easily adjusted to meet your child's needs and dysfunctions. Pelvic belt, headrest, height adjustable seat and armrests further enhance comfort and safety of the user.

It is made of user-friendly materials. Available in **5** sizes.



### IT IS RECOMMENDED IN THE FOLLOWING CASES:

- Cerebral palsy – various types and forms.
- Brain and spinal cord injuries with lower limbs paralysis and weakness, balance and motor coordination disorders with different etiology.
- Demyelinating disease (e.g. multiple sclerosis – SM).
- Meningocele.
- Muscles diseases with weakness and paralysis (dystrophy, myasthenia gravis, myositis).
- Genetic disorders.
- Degeneration and neurological diseases.
- Other condition connected with muscles weakness, paralysis and musculoskeletal disorders.
- After head and spinal cord injuries and spine trauma with paralysis and paresis.

### STANDARD EQUIPMENT:

Angle adjustable backrest | Adjustable depth of the seat | Medical Upholstery | Pelvic belt | Adjustable armrests



## ADVANTAGES:

**Stable structure**

The product ensures adequate stability and safety of a Patient.

**Eco-friendly materials**

We use only human and environment friendly materials.

**Tilt-in-space**

It is possible to change the tilt angle of the seat.

**Easy to use**

Intuitive way of adjustment and easy assembly.

**Hold and pull©**

Our harnesses and pelvic belts are equipped with Hold and pull system that provides fast and comfortable support your child in the device.

**24 months warranty**

We offer 24 months warranty on the product. We also provide after sales service.



Adjustment of a seat depth and backrest tilt angle available as a standard equipment.



Tilt-in-space function is available thanks to the legs height adjustment.



HexagonPro© Headrest has been designed to create the most comfort for user by supporting and protecting their head.



Medical upholstery used in this product is vapor permeable, waterproof and easy to disinfect.

## COLOURS

## MEDICAL UPHOLSTERY MEDICUS™



WATERPROOF



EASY TO CLEAN



STEAM PERMEABLE



ECOLOGICAL



FIRE RETARDANT



MED Sea  
Standard



MED Blue  
Optional\*



MED Skin  
Optional\*

## MEDICAL UPHOLSTERY EXTREME-MED™



BPA FREE

PHthalate FREE

300,000  
MARTINDALE  
RUBS

ABRASION



DURABLE



ECOLOGICAL



CRIB7

ANTIMICROBIAL  
& ANTIFUNGAL

STAIN RESISTANT



Citrus Green  
Optional\*



Candy  
Optional\*



Pillarbox  
Optional\*

Optional\* colors available for special request. Lead time ca. 14 days

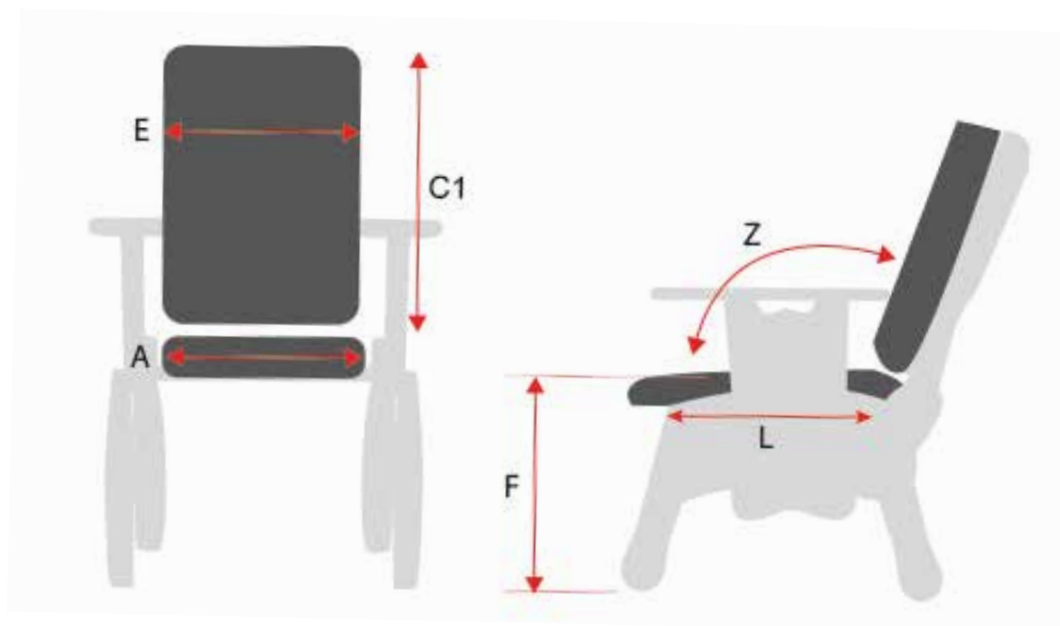
*Our passion - your independence!*

Colors may differ from the presented images. The manufacturer reserves the right to make changes in products without updating the directory.  
Current offer on our website.

Akces-MED sp. z o.o.  
Jasionka 955B  
36-002 JASIONKA

export@akces-med.com  
+48 (17) 864 04 77  
www.akces-med.com

Version: 14/05/2018



Symbol	Measurement	Unit	Size 1	Size 2	Size 3	Size 4	Size 5**
A/E	Seat/ Backrest Width	[cm]	23	26,5	31	35	42
L	Seat depth	[cm]	18-23,5	22,5-30	24,5-34,5	31-41	36-47
C	Backrest height	[cm]	26	31	35	40	44
F	Floor to seat height	[cm]	19,5-26	24-31.5	27,5-36,5	33-44	36-49
Z	Backrest tilt angle	[°]	96-116	95-116	95-116	95-116	95-116
Max. user weight		[kg]	20	30	40	55	70
Indicative user height*		[cm]	90-100	90-115	100-130	110-115	120-160

Product dimensions & weight							
Width	[cm]	34	37	42	47	55	
Length	[cm]	37.5	41	53	58	66	
Height	[cm]	48	56	63	75	81	
Weight	[kg]	8	10	12	15	14	

**Note:** Measurement tolerance range is +/- 1 cm / ½ in.

\*The value cannot be main factor for selection of product size.

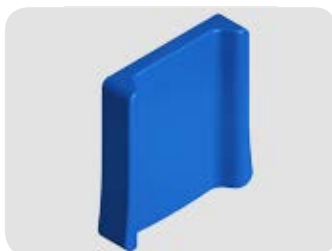
\*\*For special request.



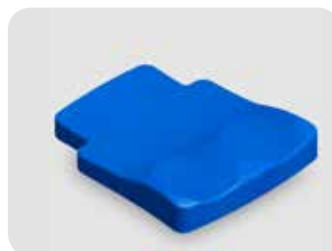
## Comfort

**Tray**

Easy to assemble / disassemble. It is also child-proof, which means that your child cannot release the tray. Light and easy to clean. Perfect for playing and feeding. Provides also extra trunk support

**Shaped backrest cushion**

Backrest cushion with an ergonomic shape which has a significant impact on the patient's comfort while using the chair.

**Shaped seat cushion**

Seat cushion with anatomic cuts for thighs which provides a comfortable position while using the chair.

**Table JOY™**

School table that can be used with Jordi or any other school chair. Height adjustable, easy to clean and very durable. JOY™ table is available in various options.

**JOY™ table is available in following options:**

**436/1 – Cut-Out Board**

**436/2 – Plain Board**

**436/3 – Tilt Cut-Out Board**

**436/4 – Tilt Plain-Board**

**Feet Cover**

Comfortable pad, protecting the feet from direct contact with the footrest



## Positioning

**Hexagon Pro™ headrest**

Headrest adjustable in three dimensions. Due to its special design, it can be precisely adjusted to the user's needs.

**Knee abductor**

Adjustable and easy to remove. Prevents in legs crossing. Made of soft foam, upholstered with easy to wash medical fabric. Collision with abduction block.

**Narrowing pads (2 pcs) 3cm**

Made of elastic foam trimmed with medical upholstery. They are attached to the armrests from both sides. Thanks to them you can narrow the seat by 6cm.

**Narrowing pads (2 pcs) 5cm**

Made of elastic foam trimmed with medical upholstery. They are attached to the armrests from both sides. Thanks to them you can narrow the seat by 10cm.

**Laterall supports**

Width adjustable supports that provide additional side stabilization of the trunk. You can mount them at different height.

**Abduction block**

Shaped abduction block provides with correct seating position. Can be easily removed. Collision with: Knee abductor. Restriction: tray height

**4 points safety vest**

It acts as a trunk stabilizer. Wide range of adjustment helps to fit it to the user's needs. Equipped with HOLD&PULL® system.

**Abduction-adduction thigh belts**

They can be easily attached to the chair by use of Velcro straps. Their use allows to keep legs in abduction and also prevents the patient from slipping off of the device.



**Safety****Skis / rocking chair**

Helps to stabilise and prevent tip over and chair movements. Once revert the skis upside down Kidoo become a rocking chair.

**Build-up****Footrest**

Increases the comfort of using the chair and provides foot support. You can adjust the height.

**Mobile Base**

You can easily transfer Kidoo from one place to another.

**Mobile base with footrest with 2 points foot stabilizing belts**

Provides foot support and offers mobility.

**Footrest with 2 points foot stabilizing belts**

Can be attached to the Kidoo chair or it can act as an independent footrest. Height and angle adjustable making it easy to tailor to the needs of your child. Collision with: skis, mobil base, mobil base with footrest,

**Mobile base with footrest**

Mobile base with footrest for easy movement.