

ABCO Adult Change Table




ABCO Health Care's Adult Change table is designed and manufactured in Australia for use in schools, disabled accommodation, sports & special needs facilities and in other locations where clients need help with toileting or changing. With a full-sized table top it is long enough and wide enough for most patients. With a minimum height of less than 500mm clients can transfer to the change table independently, from a wheel chair or with a hoist [there is 155mm clearance below the frame for hoist clearance]. Carers can adjust the working height to the best level for them [up to 950mm] using the foot control. Accessory generously upholstered cot sides and cot ends can also be fitted for extra client security. The four individually braked wheels make moving the table easy without compromising stability in. One wheel has a directional locking brake to help control the table when moving. There is no Australian GST payable on Change Tables.

Shown with optional upholstered Cot Side

Product Code: **CA00.07021.E0000F.00C10**

Standard Configuration	Factory Fitted Options
710mm wide Upholstery	800mm
200kg Lifting Capacity	350kg [halves lifting time]
Wedgewood upholstery	Other vinyl colours
Foot control	Hand control
	Add Battery
	Add Upholstered Cot Side [per side]
	Add Upholstered Cot End [per end]
	Strap Bar [included with Cot Sides][per side]
	Patient Strap [Also order Cot side[s] and/or Strap bar[s]]

Detailed Standard Specification		
Couch Dimensions	Overall Length	1930mm
	Width	710mm
Upholstery	Fabric	Vinyl – Medical grade heavy duty vinyl fabric. Formulated to inhibit hardening and breakdown of the PVC surface due to perspiration and oils. Light fast, fire retardant, sanitised and stain resistant.
	Foam	38mm high density foam rolled over the edges of the strong plywood base covered with 12mm of soft foam for enhanced comfort while providing a full working width and extended life.
	Colours	Wedgewood. Other colours on request
Frame	Materials	Heavy steel plate and tube
	Colour	Light Grey
	Wheel System	125mm braked castors providing 155mm clearance below the frame for patient hoist access.
Height adjustment	Range	450 – 950mm
	Lift mechanism	Two leg swing arm
	Lift operation	Linak brand European electric actuator
	Control	Foot control
Electrical Requirements		Medically approved to AS 3200 / EN60601 - 1
	Supply	230vAC $\pm 10^\circ$, 50hz max 1.5amp
	Operation	24vDC max 70VA
Lifting Capacity		200kg
Weight		75kg
Shipping Dimensions		197 x 80 x 60cm / 90kg. Shrink-wrapped on a wooden pallet with bubble wrap
ARTG No		123753
Warranty		5 years, parts only

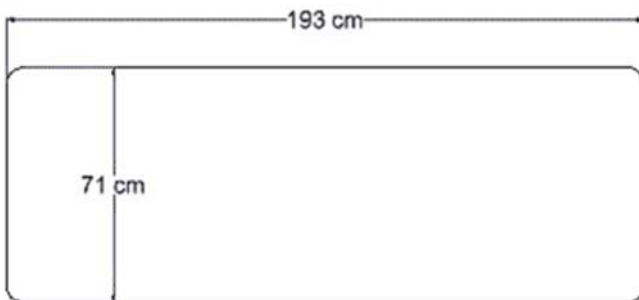

Height Range 45 - 95cm

ABCO Health Care is a specialist manufacturer of electrically operated couches, tables and chairs for health care applications.

This brochure shows only part of our extensive range. Please visit our web site at www.abcohc.com.au

ABCO Health Care has a policy of continuous improvement. Specifications and prices are subject to change without notice.

ABCO Health Care Pty Ltd

A.B.N. 71 073 519 356

Australian Designed, Made and Owned.

3 – 7 Discovery Road,
Dandenong South VIC 3175 Australia

Phone +61 03 9768 3756

E-mail: sales@abcohc.com.au